

# Nonin CO-Pilot™

## Wireless Handheld Multi-Parameter System (H500)

---

### Quick Start Guide






Refer to the CO-Pilot Wireless Handheld Multi-Parameter System (H500) Operator's Manual for complete instructions for use to fully understand the device, its operation, and the potential risks to user and patient.

 Operator's Manual included on provided USB flash drive

## Table of Contents

01	System Components
02	Display Icons and Indicators
03	Initial System Set Up
04	Taking a Measurement
05	Placing Disposable rSO <sub>2</sub> Sensor
06	Placing Finger Clip Sensor
07	Determining Pulse Signal Quality
08	Checking Signal Processor Battery Status
09	Charging the System
10	Long Term Storage Mode
11	Basic Troubleshooting
12	COHb Quick Reference
13	COHb Exposure Treatment Decision Process

## 01 System Components

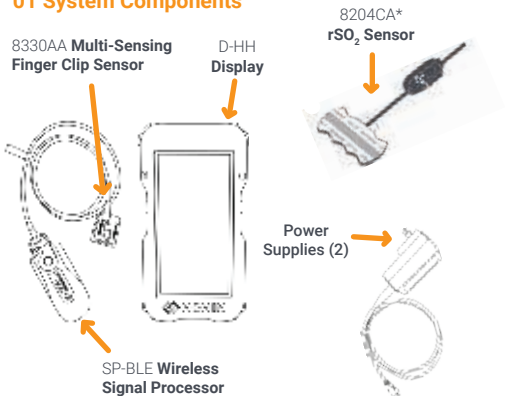

\*rSO<sub>2</sub> sensors not included in standard configuration  
\*\*Not depicted: Data Download Cable and USB Flash Drive

## 02 Display Icons & Indicators

### Multi-Sensing Display Screen



Display Battery Indicator

Pulse Signal Quality Indicator

Pulse Wave Form/ PPG

Display Charging Port

### rSO<sub>2</sub> Display Screen



rSO<sub>2</sub> Histograms  
Trending readings of the rSO<sub>2</sub> value over time

On/Off, Sleep Button


Session ID

Signal Processor Serial Number

Event Marker Button

Signal Processor Battery Life Indicator

### Signal Processor (SP-BLE)

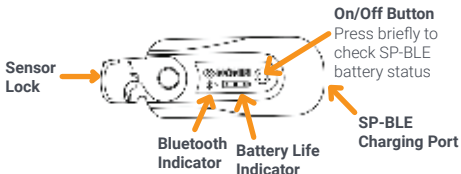

On/Off Button  
Press briefly to check SP-BLE battery status

Sensor Lock

Bluetooth Indicator

Battery Life Indicator

SP-BLE Charging Port

## 03 Initial System Setup

- 1 Charge display and signal processor per instructions in Section 9
- 2 Press and hold On/Off button on signal processor until green LEDs blink on and blue Bluetooth LED begins to blink
- 3 **Press and hold** On/Off, Sleep button on display until screen turns on
- 4 Connect sensor to signal processor and engage sensor lock (Figure 1)
- 5 System is now ready to take a measurement or to be put into Sleep Mode for future use

Figure 1



## 04 Taking a Measurement with CO-Pilot

- 1 Activate the display from sleep mode by briefly pressing the On/Off, Sleep button on the right side of the device (Figure 2)
- 2 When the display is activated, it will automatically be connected to the signal processor and the sensor will turn on
- 3 If the signal processor is not automatically activated, press and hold down the On/Off button on the signal processor
- 4 Place the sensor on the patient according to sensor instruction
- 5 When operation is complete, briefly press the On/Off, Sleep button on the display to end the session and put device into Sleep Mode; signal processor will enter Sleep Mode automatically



Press On/Off,  
Sleep button  
to activate  
display  
Figure 2

## 05 Placing Disposable rSO<sub>2</sub> Sensor

- 1 Connect INT-100 to signal processor as shown in Figure 1
- 2 Connect 8204CA sensor to INT-100 and engage lock per Figure 3
- 3 Remove the protective backing from the sensor pad and gently, but firmly place the sensor above the eyebrow per Figure 4
- 4 Ensure sensor surface adheres to the skin to prevent light from traveling between emitting or receiving elements or ambient light from entering



Figure 3

### NOTE

An improperly placed sensor may result in inaccurate readings. If sensor does not adhere fully, sensor may result in inaccurate reading



Figure 4

## 06 Placing Finger Clip Sensor

- 1 Connect sensor to signal processor as shown in Figure 1
- 2 Select proper sensor size and use a clean, dry finger. Heavy lotions, sunscreen, dirt, oils and nail polish may alter accuracy
- 3 Gently squeeze open the sensor and insert the finger all the way into the sensor

## 07 Determining Pulse Signal Quality

Pulse signal quality can be determined by the color of the pulse waveform as shown.



Green = Good



Yellow = Marginal



Red = Insufficient

## 08 Checking Signal Processor Battery Status


To check the signal processor battery status, quickly press the On/Off button on the signal processor. Green battery indicator lights will light indicating current battery status. Signal processor battery status can also be found on the display screen as shown in Section 02.

### Signal Processor Battery Status Matrix

Battery Indicator	Battery Level
Three green lights lit	High (60-100%)
Two green lights lit	Medium (20-59%)
One green light lit	Low (5-19%), charge SP-BLE
One green light blinking	Critical (0-4%), charge SP-BLE
No Lights	Depleted or off

## 09 Charging the System

- 1 Charge the display and signal processor using the provided Nonin power supplies **ONLY\***
- 2 Insert power supply into respective ports shown below to charge.

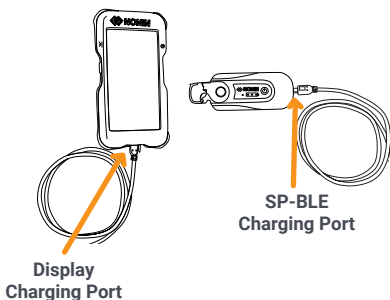

## 10 Long Term Storage Mode

- 1 To place the display into long term storage mode, **press and hold** the On/Off, Sleep button on the display (Figure 2) until screen prompt appears
- 2 Choose "Power Off" on display screen. Screen will go black, display is now in long term storage mode
- 3 To place the signal processor into long term storage mode, press and hold the on/off button on the signal processor until all lights flash simultaneously. Signal processor is now in long term storage mode



Figure 5

## 11 Basic Troubleshooting



### Sensor Not Connected to Signal Processor

Check that sensor is plugged into signal processor and sensor lock is engaged



### Sensor Fault

Sensor may be inoperable, try another sensor



### Inadequate Signal

Ensure sensor is properly applied to patient



### Cannot Connect to Signal Processor

Check that signal processor is on, battery is charged, and in range of display



### Error Code

Refer to Operator's Manual

## 12 COHb Quick Reference

Carbon monoxide (CO) is an odorless, colorless gas that can remain undetectable until exposure results in injury or death. CO poisoning can lead to a variety of clinical effects that resemble many neurological and cardiovascular conditions, leading to potential misdiagnosis.

### Common Symptoms of CO Poisoning

- 1 Headache
- 2 Dizziness
- 3 Fatigue and Weakness
- 4 Loss of Consciousness
- 5 Nausea

### Important Symptom Considerations

- A patient exposed to high concentrations for a short time may be less symptomatic than a patient who reaches the same COHb level after a prolonged exposure.
- Chronic poisoning is generally associated with the less severe symptoms.
- Some patients presenting with a carboxyhemoglobin level of 20% may be remarkably symptomatic, while others experiencing the same level of COHb% may exhibit only mild, equivocal symptoms.
- Low-level exposure can exacerbate angina and chronic obstructive pulmonary disease.
- Coronary artery disease patients are at risk for ischemia and MI even at low levels of CO.
- Severity of clinical condition is related to CO concentration, the duration of exposure, and the prevailing clinical disposition of the patient.
- Smokers present with higher levels of COHb than do non-smokers. The COHb level in non-smokers is approximately 1-2%. In smokers, a COHb of 3-8% would be considered normal, depending on severity of smoking.



## Treatment

Carbon Monoxide: Half Life Elimination from Blood	
Room Air	240-360 minutes
Oxygen (100%)	80 minutes
Hyperbaric Oxygen (HBO)	22 minutes

## 13 COHB Exposure Treatment Decision Process

If CO exposure is NOT suspected:

- Noninvasive COHb values should NOT be used

If CO exposure IS suspected:

- Patient has signs and symptoms, transportation to the hospital is recommended even if the COHb is <10%
- Patient has COHb  $\geq$ 10%, transportation to the hospital is recommended, even if asymptomatic



## Nonin Medical, Inc.

13700 1st Avenue North  
Plymouth, MN • 55441-5443 • USA  
**Tel:** +1.763.553.9968 1.800.365.8874  
**Fax:** +1.763.577.5521  
**Email:** info@nonin.com

## Nonin Medical B.V.

Doctor Paul Janssenweg 150  
5026 RH Tilburg • Netherlands  
**Tel:** +31 (0)13-45 87 130  
**Email:** infointl@nonin.com



<https://www.nonin.com/supr/orioz-pilot>

**ENG** - Translations of this quick start guide can be found using this QR code or on the USB for the CO-Pilot™ operator's manual.

**FRE** - Vous pouvez trouver les traductions de ce guide de démarrage rapide à l'aide de ce code QR ou sur la clé USB pour le manuel de l'opérateur du CO-Pilot™.

**GER** - Übersetzungen dieser Schnellstartanleitung finden Sie über diesen QR-Code oder auf dem USB-Stick mit der CO-Pilot™-Bedienungsanleitung.

**ITL** - Le traduzioni di questa guida rapida possono essere trovate utilizzando questo codice QR o sull'USB per il manuale dell'operatore CO-Pilot™.

**SPA** - Las traducciones de esta guía de inicio rápido se pueden encontrar usando este código QR o en el USB del Manual del operador de CO-Pilot™.

**POR** - É possível obter as traduções deste guia de início rápido utilizando este código QR ou na unidade USB com o manual do operador CO-Pilot™.

**DUT** - Vertalingen van deze beknopte gebruiksaanwijzing zijn te vinden via deze QR-code of op de USB-stick voor de CO-Pilot™-gebruikershandleiding.

**GRK** - Μπορείτε να βρείτε μεταφράσεις αυτού του οδηγού γρήγορης εκκίνησης χρησιμοποιώντας αυτόν τον κωδικό QR ή στο USB για το εγχειρίδιο χειριστή CO-Pilot™.

**DAN** - Oversættelser af denne startvejledning kan findes ved hjælp af denne QR-kode eller på USB'en til CO-Pilot™-brugervejledningen.

**SWE** - Översättningar av den här snabbstartsguiden hittar du via den här QR-koden eller på USB-enheten för CO-Pilot™ användarhandbok.

**FIN** - Tämän pika-aloitusoppaan käännökset löytyvät tästä QR-koodista tai CO-Pilot™-käyttöoppaan USB-muistitikulta.

**POL** - Tłumaczenie niniejszej skróconej instrukcji obsługi można znaleźć, skanując kod QR lub na nośniku USB instrukcji obsługi urządzenia CO-Pilot™.

**NOR** - Oversettelser av denne hurtigstartveiledningen finner du ved hjelp av denne QR-koden eller på USB-en for bruksanvisningen til CO-Pilot™.