

INSTALLATION GUIDE

This document explains how to download, install, and register the nVISION® Data Management software on a Windows PC.

System requirements necessary to run nVISION software:

- nVISION Data Management Software is compatible with Windows XP SP3, Windows 7, Windows 8, Windows 10 (32- and 64-bit), Windows 11 (64-bit)
- At least one serial communications port
- 256MB memory space
- 20MB free space on the hard drive
- Administrator's rights
- Internet access
- Registration key code

Step 1: Request access to download nVISION at download.nonin.com

1. Open the software download email from Nonin Medical.
2. Click **"DOWNLOAD your software now."** This link from the order shipment email will automatically take you to the download website at download.nonin.com.

Figure 1: Email confirmation showing the nVISION download link



Hello Emma,

Your digital download request has shipped. (Order#: 6810 placed on 07/02/2026)

[Click here to DOWNLOAD your software now.](#)

You may also use the preceding link to review your order details or print an invoice

nVISION® 6.5.3 Data Management Software Upgrade - Download and Installation Instructions

1. Select - Click here to download your software now.
2. Select - Click here to download the software you've purchased.
3. Click - nVisionVer 6_5_3.exe
4. Click - Run, you may see a Windows warning at the bottom of your screen, Click - Actions, Click - Run anyway
5. Click - Extract
6. nVISION - InstallShield Wizard will run

Thanks again for shopping at [Download.nonin.com!](https://download.nonin.com)

Visit us anytime at <http://hdaby-amqcp.volusion.store/>

Step 2: Access your download by navigating to the Nonin portal

1. Enter your email address and password.
2. Click **Login**.

Figure 2: Nonin login page

Step 3: Locate your order

1. Navigate to **My Account**.
2. Select **My Orders - > Review order/track packages**.
3. Locate the nVISION order.

Figure 3: Order History page

Step 4: Download the installation file

1. Open the order details. Click **View Order or Download File(s) >>>**
2. Under **Your Downloads**, click here to download the files you've purchased.
 - o **nVISIONVer6_5_3_0.exe**

Figure 4: Order Details page with download link

The screenshot shows the NONIN nVISION user interface. At the top, there's a header with the NONIN logo and 'nVISION My Cart (0)'. Below this is a 'MY ACCOUNT' section with a 'Return to My Account' link. The main area is titled 'Order History' and includes a search bar with the text 'Don't see what you're looking for? Choose the time period you wish to view' and a dropdown menu set to 'Orders placed in the past 1 month'. A 'GO' button is below the dropdown. Below the search bar is a table titled 'My Orders'. The table has columns for 'Order date', 'Items Ordered', 'Price', and 'Status'. The first row shows an order from 07/02/2026 for '1 of : nVISION V6.5.3 Upgrade' at a price of '\$0.00' with a status of 'Shipped'. To the right of the table is a 'View Order' button and a link 'Download File(s) >>>'. Below the table is a section titled 'Order Details' with a link 'Click here to download the files you've purchased.' and a checkbox 'Need to print an invoice?'. The order details section shows 'Order#: 6810' and 'Order Placed: 07/02/2026 at 08:43AM PST'. Below this is a table titled 'Order Status: Shipped' with columns for 'Ship To', 'Items Ordered', 'Price', and 'Status'. The 'Ship To' column shows the address: 'Al-Shifa Clinic, Noor Mujir, 1401 Gardena Avenue, Fridley, MN 55432, United States, 6123149344'. The 'Items Ordered' column shows '1 of : nVISION V6.5.3 Upgrade'. The 'Price' column shows '\$0.00' and the 'Status' column shows 'Shipped'. Below the table is a section titled 'Payment Information' with a 'Bill To' address: 'Al-Shifa Clinic, Noor Mujir, 1401 Gardena Avenue, Fridley, MN 55432, United States, 6123149344'. The 'Payment' column shows 'Payment:'. The 'Order Summary' section shows 'Product Subtotal: \$0.00', 'Tax: \$0.00', and 'Shipping & Handling: \$0.00'.

Step 5: Run the download installer

1. Open Chrome Downloads. Click **nVisionVer6_5_3_0.exe**. If there is a security prompt, click ... More Actions next to the file name and select Keep.
2. Click the **down arrow** next to delete and select **Keep anyway**.
3. Locate and open the **nVISIONVER6_5_3_0.exe**. file.

Figure 5: Downloads page showing the installer file Details page with download link


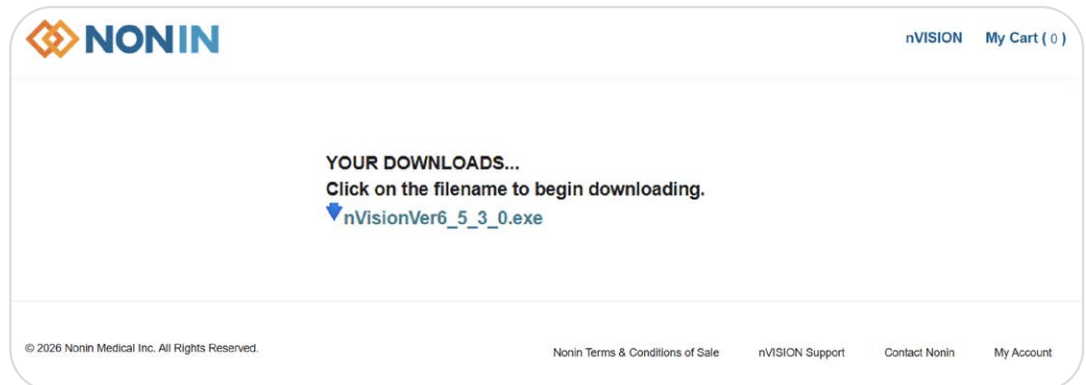


Figure 6: Windows downloads pop-up showing the link for the installer file

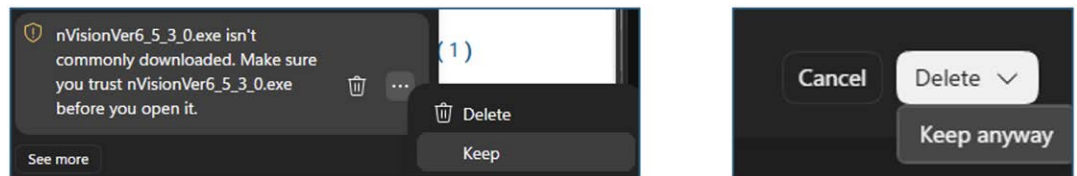

Step 6: Launch InstallShield Wizard and accept the License Agreement

1. After extraction, when the InstallShield Wizard starts, Click **Next**.
2. Review the license agreement. Select: **I accept the terms of the license agreement**.
3. Click **Next**.

Figure 7: InstallShield Welcome screen


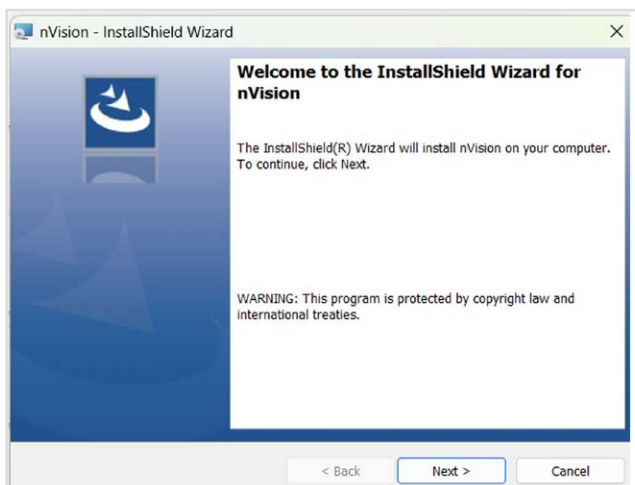


Figure 8: License Agreement screen

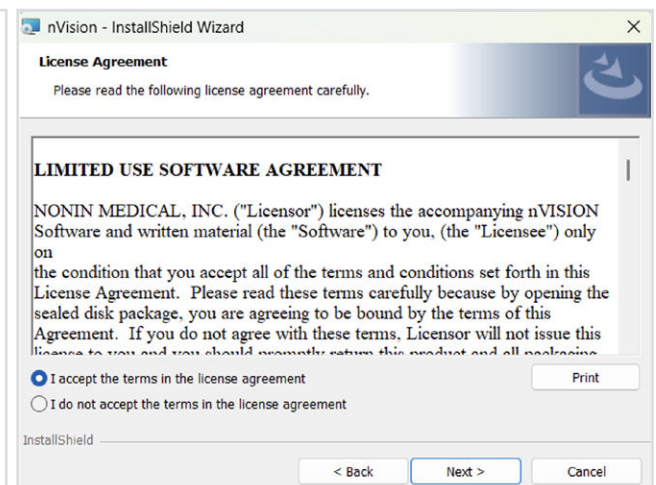

Step 7: Start the installation after selecting the location

1. Review the destination folder: C:\Program Files (x86) \Nonin Medical, Inc.\nVISION\
2. Click **Next**.
3. Review installation settings and click **Install**.

Figure 9: Destination Folder screen


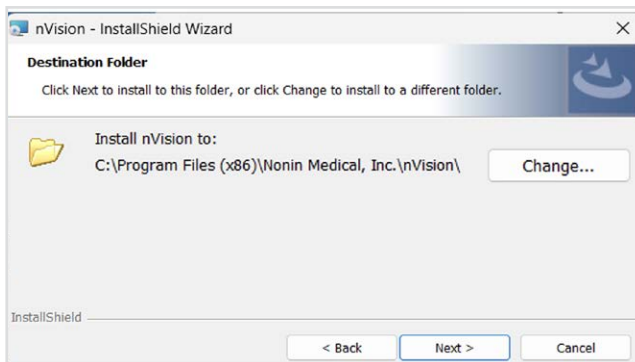




Figure 10: Ready to Install screen

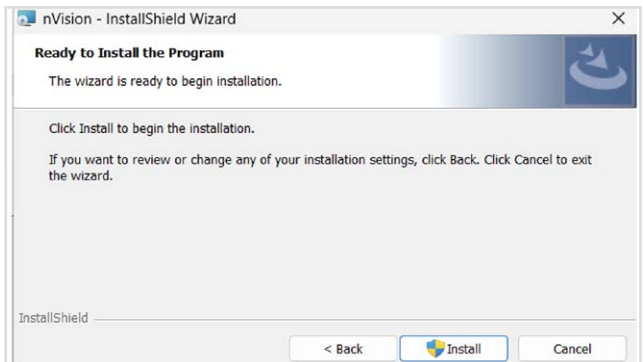

Step 8: Approve user account control

1. When prompted by Windows, click **Yes**.

Figure 11: User Account Control prompt

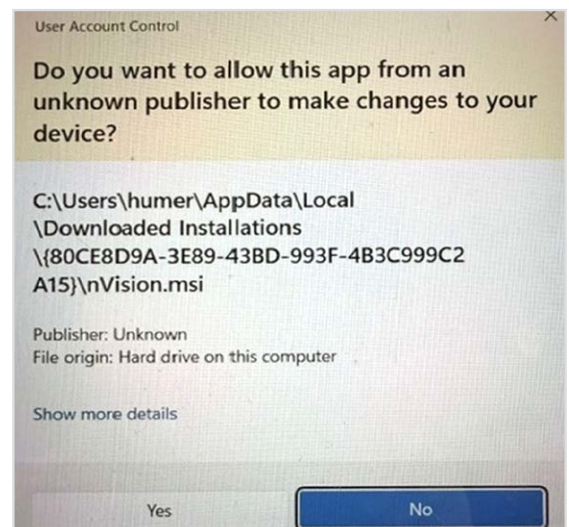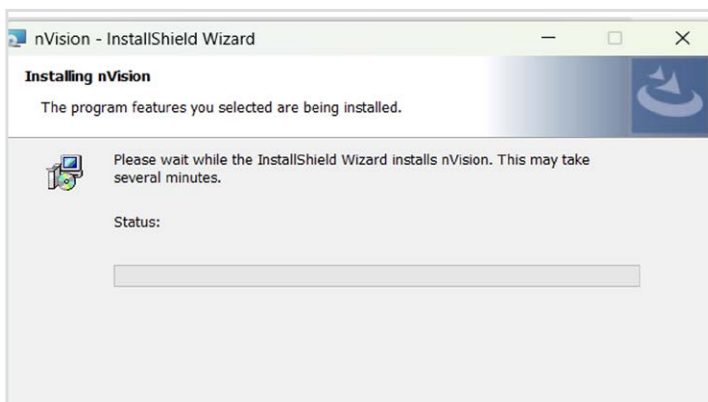

Note: Windows SmartScreen may display a warning because the installer has not yet been recognized by Microsoft's reputation-based protection system. Verify that the file was downloaded from the official Nonin portal before proceeding.


Step 9: Complete installation and register your software

1. Wait for installation to finish.
2. Select **Launch nVISION** if desired.
3. Click **Finish**.
4. Enter: Registration Information
 - a. User Name
 - b. Company
 - c. Key Code (must include dashes)
5. Click **OK**. The home screen is displayed

Figure 12: InstallShield Wizard Screen



Figure 13: nVISION Registration Window

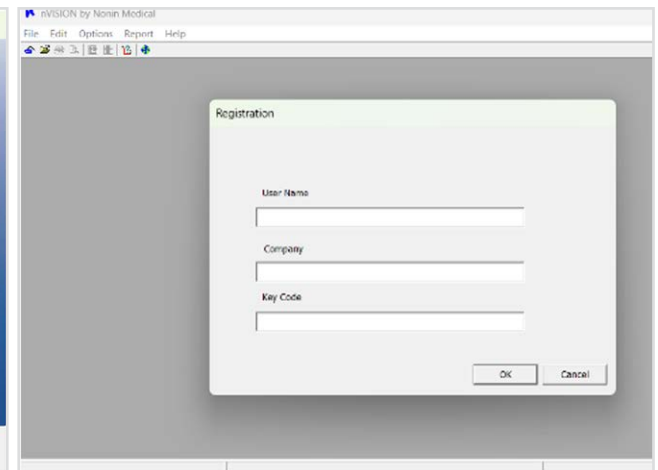

Step 10: Verify successful installation

1. Launch nVISION.
2. Confirm the application opens without errors.
3. Verify that the registration information is accepted.
4. Confirm that nVISION software version 6.5.3 is installed.

Troubleshooting Tips

Common Issues

1. **Issue:** Download link does not open
Resolution: Verify login credentials and internet connection.
2. **Issue:** SmartScreen blocks installation
Resolution: Click **...More Actions** next to the download file name and select **Keep**.
Click the down arrow next to delete and select **Keep anyway**.
3. **Issue:** Registration key rejected
Resolution: Verify User Name, Company, and Key Code entries.

For troubleshooting assistance or any questions regarding installation, please contact **Technical Support Service** at technicalservice@nonin.com.

Nonin Medical, Inc.

13700 1st Avenue North, Plymouth MN • 55441-5443 • USA

Tel: +1.763.553.9968 OR +1.800.365.8874

Email: info@nonin.com



nonin.com