



Model 7500 Quick Start-up Guide

SET USER-DEFINED DEFAULTS

SET PATIENT SECURITY MODE ON

DATA COLLECTION & DOWNLOAD

TO CLEAR MEMORY

1. Turn unit on – Review factory settings below
2. Step through ALL settings using LIMITS button
3. Change desired settings using or buttons
4. MUST push LIMITS button after each setting to confirm & move on
5. Scroll to end of settings to 2nd symbol (=minutes) and press LIMITS button
6. Unit will automatically review new settings
7. TO SAVE SETTINGS – press and hold ALARM SILENCE button while briefly pressing LIMITS button – “dEF On” will flash

Return to Factory Defaults

Note: User-Defined Default values are lost when returned to Factory Defaults.

1. Turn unit on – “dEF On” will appear
2. Press and hold ALARM SILENCE button and briefly press button – “dEF OFF” will flash
3. Press POWER button until unit turns OFF, then turn back ON

Factory Default Settings

The values at right come pre-programmed in the Model 7500. The Patient Security Mode is OFF, unless turned on.

KEY	Alarm Limit	Factory Setting	Adjustment Options	Increment
	SpO ₂ Low Alarm Limit	85%	Off, 50-95	1%
	Pulse Rate High Alarm Limit	200 BPM	Off, 75-275	5 BPM
	Pulse Rate Low Alarm Limit	50 BPM	Off, 30-110	5 BPM
	SpO ₂ High Alarm Limit	Off	Off, 80-100	1%
	Alarm Volume	High	Off, Low, High	N/A
	Oximetry Averaging	No	Yes, No	N/A

1. Turn unit on – Review factory settings below
2. Step through ALL settings using LIMITS button
3. Change desired settings using or buttons
4. MUST push LIMITS button after each setting
5. Scroll to end of settings to 2nd symbol (=minutes) and press LIMITS button
6. Unit will automatically review new settings
7. TO SAVE SETTINGS – press and hold ALARM SILENCE button while briefly pressing LIMITS button – “dEF On” will flash
8. Press POWER button until unit turns OFF
9. TO SET – Press and hold ALARM SILENCE button while turning the unit back ON
10. Screen will flash “SEC On”

Note: To set “Low O₂” alarm below 85%, must set in User-Defined Defaults first, then follow steps above.

Turn Off Patient Security Mode

1. Press POWER button until unit turns OFF
2. Together hold down ALARM SILENCE button, then LIMITS button and then POWER button
3. Units reads “SEC OFF”

Data Collection

1. Turn unit on – Recording begins automatically
- Note:** Data less than 1 minute is not saved.

Downloading Data To PC

Note: Patient Security Mode must be OFF


1. Press POWER button until unit turns OFF
2. Connect unit to PC using 7500SC download cable
3. Boot up nVISION and select FILE>New Data Capture, click 7500 oximeter option
4. On 7500 unit press and hold and POWER buttons - unit will display “PLA bAC” during download
5. Unit returns to normal operating mode after download

1. Turn unit on

Note: Patient Security Mode must be OFF

2. Press LIMITS button until “CLr no” is displayed
3. Press button until “CLr YES” appears
4. Press LIMITS button to confirm
5. “dEL no” appears, press button so “dEL YES” appears
6. Press LIMITS button to confirm
7. “dnE CLr” confirms memory has been cleared

For complete specifications, please refer to the Operator’s Manual.

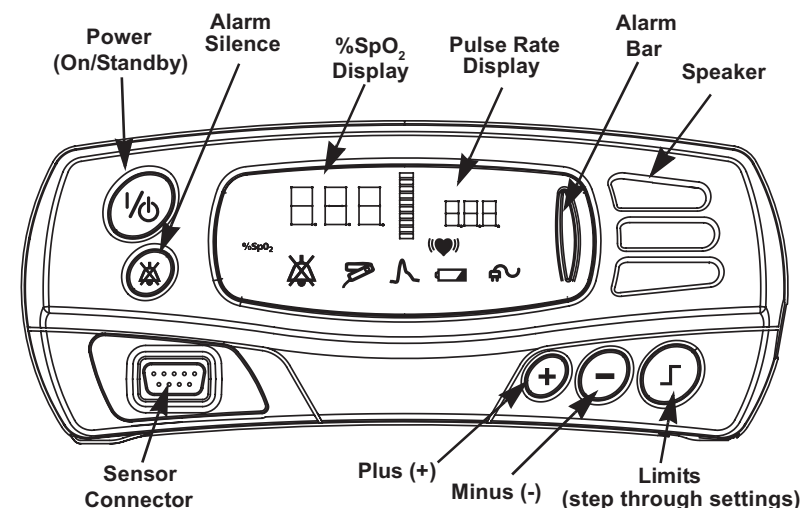



<https://www.nonin.com/support/7500/>

Model 7500

Quick Start-up Guide

EN: Translations of this quick start guide can be found using this QR code or on the USB for the Model 7500 operator's manual.

FR: Les traductions de ce guide de démarrage rapide sont accessibles à l'aide de ce code QR ou sur la clé USB contenant le manuel de l'opérateur du modèle 7500.

DE: Übersetzungen dieser Schnellstartanleitung finden Sie über diesen QR-Code oder auf dem USB für die Bedienungsanleitung des Modells 7500.

IT: Le traduzioni di questa guida rapida sono reperibili utilizzando questo codice QR o sull'USB per il Manuale dell'operatore del modello 7500.

ES: El CD del manual del operador del modelo 7500 contiene traducciones de esta guía.

PT: As traduções deste guia de início rápido podem ser encontradas ao utilizar este código QR ou no USB do manual de operação do Modelo 7500.

NL: Vertalingen van deze snelstartgids zijn te vinden met behulp van deze QR-code of op de USB voor de gebruikershandleiding van Model 7500.

EL: Το εγχειρίδιο χειριστή για το WristOx2®, Μοντέλο 3150, βρίσκεται στη μονάδα USB που παρέχεται ή online χρησιμοποιώντας τον παραπάνω κωδικό QR.

DA: Oversættelser af denne startvejledning kan findes ved hjælp af denne QR-kode eller på USB-drevet med brugervejledningen til Model 7500.

SV: Översättningar av denna snabbstartsguide kan nås med QR-koden eller på USB-enheten med bruksanvisningen för modell 7500.

FI: Tämän pikaoppaan käännökset löytyvät tämän QR-koodin avulla tai USB-muistitikulta Model 7500 -käyttöoppaan osalta.

PO: Tłumaczenie tej skróconej instrukcji obsługi można znaleźć za pomocą tego kodu QR lub na napędzie USB do instrukcji operatora modelu 7500.

NO: Oversettelser av denne hurtigstartveiledningen finner du ved hjelp av denne QR-koden eller på USB-en for brukerhåndboken Model 7500.

